

The logo for Cenefom, featuring the brand name in a blue, sans-serif font. It is positioned inside a large, stylized graphic of an eye, which is composed of overlapping circles and lines in shades of orange, yellow, and blue. The background is a light gray with a subtle pattern of dashed white circles.

OPHTHALMIC Series

Cenefom provides a full range of PVA ophthalmic products are designed to fluid control and safe cleaning of microsurgical instruments during of ophthalmic procedures.



- Designed for ophthalmic surgeries
- Lint and fiber-free
- Fast wicking and high fluid retention
- Sterilized by gamma radiation
- Packaged in cleanroom facilities

www.cenefom.com

OPHTHALMIC SERIES

PVA EYE SPEAR

Cenefom PVA eye spears are constructed of soft, absorbent and highly retentive polyvinyl alcohol foam. The porosity is controlled for a variety of applications including LASIK and refractive surgeries. Cenefom PVA eye spears, gentle on the eyeball and cornea, are designed to quickly clean the surface of the eye and absorb diffuse fluids during ophthalmic procedures.

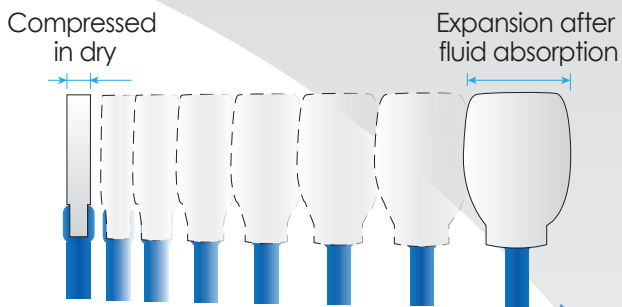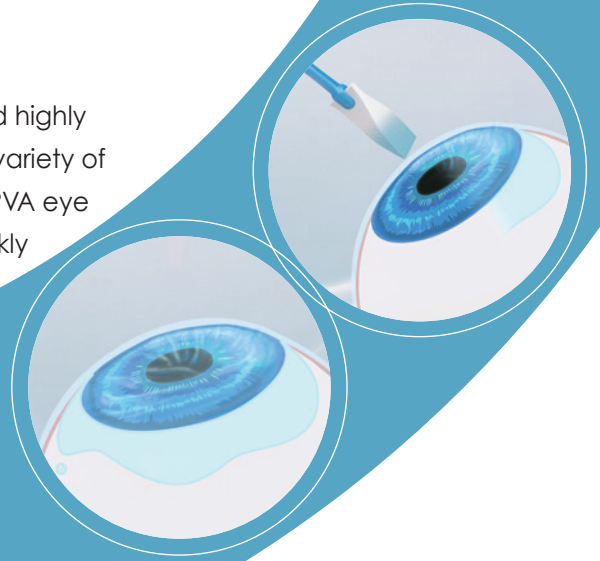


● **Fast wicking and high fluid retention**

● **Lint-free, fiber-free and particulate-free**

● **Excellent tensile strength**

● **Well bio-compatibility**

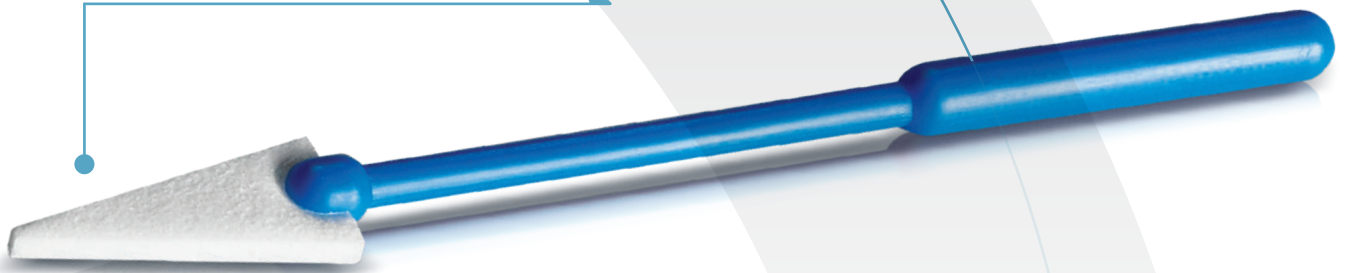

Model #	Dimensions	Packaging
OP-ES01	8 x 17 mm	5 pcs/foil bag, 2 foil bags/pouch, 20 pouches/box
OP-ES02	8 x 17 mm	5 pcs/foil bag, 80 foil bags/PE bag
OP-ES06	8 x 17 mm	5 pcs/foil bag, 1 foil bag/pouch, 40 pouches/box

OPHTHALMIC SERIES

Cenefōm eye drain comes with a fluid collection bag and an eye wick, to provide continuous control of fluid throughout the procedure. The eye wicks are made from highly absorbent, fast wicking, lint free PVA material.

EYE DRAIN & WICK



Model	Dimensions	Packaging
OP-ED01	80 cc	1 pc/pouch, 10 pouches/box
OP-ED02	400 cc	1 pc/pouch, 10 pouches/box



Model	Dimensions	Packaging
OP-EW01	3 x 220 x 2 mm	2 pcs/pouch, 20 pouches/box

Cenefōm instrument wipes are designed to provide the surgeon with superior fiber free products for fluid control and safe cleaning of microsurgical instruments, these wipes are made of the highly absorbent, fast wicking, and lint-free PVA material and cause no harm on the sensitive instruments or optical sites.

INSTRUMENT WIPE



Model	Dimensions	Packaging
OP-IW01	76.2 x 76.2 x 2 mm	1 pc/pouch, 20 pouches/box
OP-IW05	83 x 83 x 2 mm	1 pc/pouch, 20 pouches/box

Cenefōm wipe & wick includes an instrument wipe with an attached, removable eye wick.

WIPE & WICK



Model	Dimensions	Packaging
OP-WW01	76.2 x 76.2 x 2 mm	1 pc/pouch, 20 pouches/box

OPHTHALMIC SERIES

Cenefōm cleaning block is designed to provide safe cleaning of diamond knife. It's made of highly absorbent, lint-free PVA material.

CLEANING BLOCK



Model	Dimensions	Packaging
OP-CB01	45 x 15 x 20 mm	1 pc/foil bag, 1 foil bag/pouch, 10 pouches/box

After microkeratome dissection, a LASIK ring is placed on the eye before the flap is reflected. It absorbs and prevents the excess fluids and blood from encroaching into the surgical site. Cenefōm LASIK rings also minimize the risk of under-correction and the contamination of the stromal bed by fluids from fornices.

LASIK RING



Model	Dimensions	Packaging
OP-LR01	Ø10 x Ø14 x 2 mm	1 pc/pack, 1 pack/pouch, 20 pouches/box
OP-LR02	Ø12 x Ø15.4 x 2 mm	

Cenefōm LASIK eye drains allow for the simultaneous absorption of excess fluids and protection of flap and hinge.

LASIK EYE DRAIN



Model	Dimensions	Packaging
OP-LRD01	Ø10 x Ø40.6 x 2 mm	1 pc/pack, 1 pack/pouch, 20 pouches/box
OP-LRD02	Ø12 x Ø48.9 x 2 mm	

Cenefōm corneal light shields keep the cornea moist, protect the retina from potentially harmful light during ophthalmic procedures, and enhance the probability and speed of full recovery after surgery.

CORNEAL LIGHT SHIELD



Model	Dimensions	Packaging
OP-CS01	Ø7 x 2 mm	1 pc/pack, 1 pack/pouch, 20 pouches/box
OP-CS02	Ø9 x 2 mm	